


Syringe Pump Syringe Calibration

March 2019


- Power pump on
- If not connected to power press STOP to silence the 'USEBATTERY' alarm
- Press 8 to go to 'Settings'



- Press 2 for 'Syringe Settings'

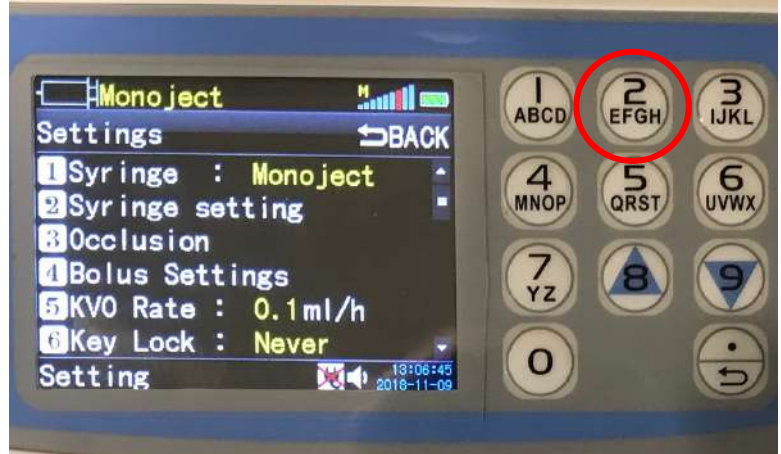

- Press 2 for 'Calibration'

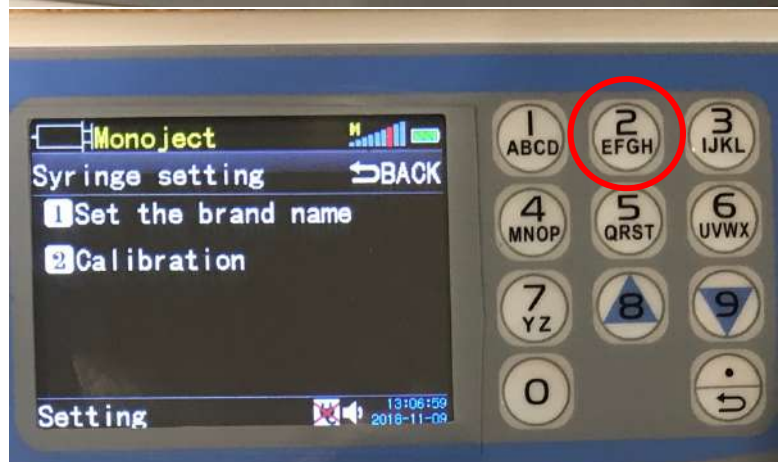

- Press 9 to scroll down to next menu




- Press brand of syringe you want to calibrate
- For this example we will be calibrating a Monoject 12 ml syringe
- Press 5 for 'Monoject'


*Nipro is no longer
programed in updated
Syringe Pumps



- Press syringe size you want to calibrate
- Press 2 for 'Set 12 ml Syringe'

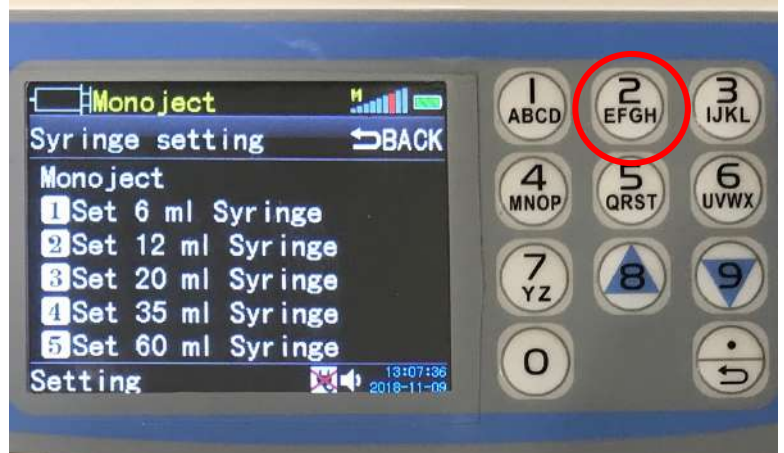

- Press 1 for 'Auto Measure Syringe'

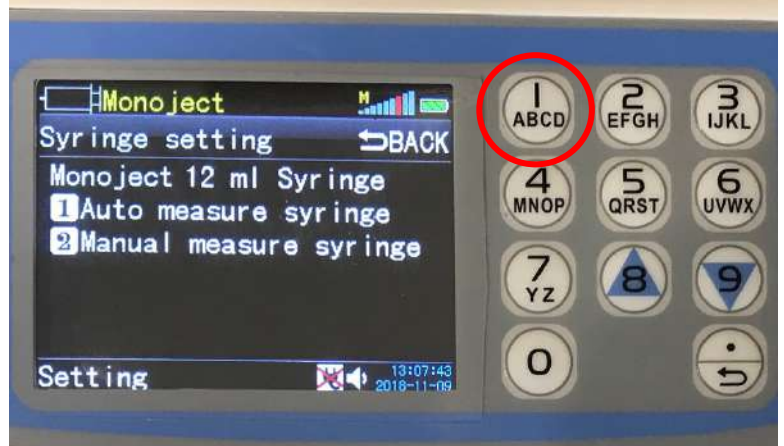

- This screen will appear




- Pull the plunger of the syringe to the 10ml graduation
- **PLEASE NOTE:** Pull the plunger to the graduation it instructs on the screen
- In this case it is a 12 ml syringe, but we pull the plunger to the 10 ml graduation

Why do you only instruct to pull to 10ml not 12ml
Need to provide an explanation.



- Install the syringe in the pump
- Once installed press 1 for 'YES'

Need to instruct the customers to
Make sure the syringe is locked
into the holder (back end)

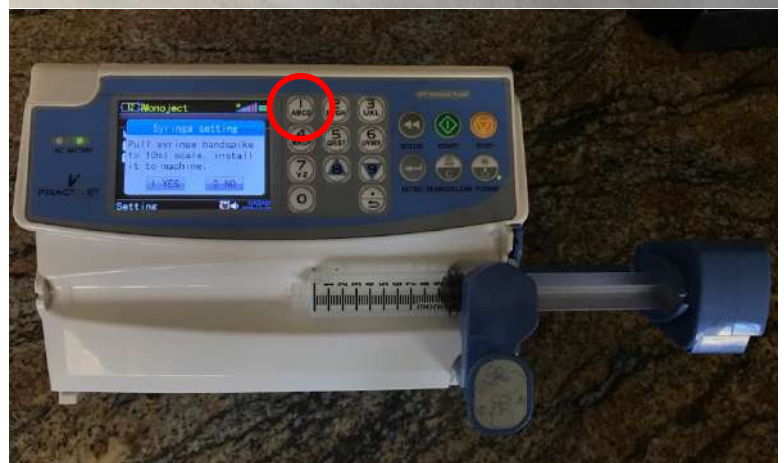

- The pump will automatically measure the diameter, barrel length, & tail length of the syringe
- This can take up to 1 minute



- Once autocalibration is done press 4 for 'Confirm'
- This syringe is now calibrated to infuse accurately
- Press the back key until you get to the main menu

