NOTE: For access to the full user manual for this product, please visit www.heska.com/productmanuals.

Vet/IV™ 2.2 Infusion Pump Overview

Fluid Choices:

- Electrolyte solutions
- Glucose-containing solutions
- Total Parenteral Nutrition (TPN)

Do not use to administer plasma or whole blood

IV Set Choices:


Hospira calibration:

- ICU Medical Brand
- Zoetis Series A
- Covetrus Series A

Baxter calibration:

- Baxter Brand
- Zoetis Series B
- Covetrus Series B

Use calibrated brand only.
Do not interchange.
Use standard sets for MACRO mode; micro sets for MICRO mode.

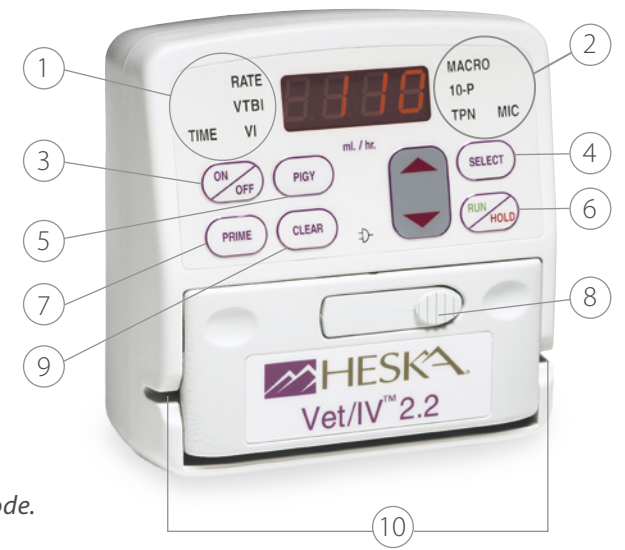

Front Panel

1. User-definable parameters:
   - Rate of infusion (RATE)
   - Volume to be infused (VTBI)
   - Time of infusion
2. Infusion mode options:
   - MACRO
   - 10–P
   - MICRO
   - TPN
3. ON/OFF button
4. Select infusion mode, pressure
5. Piggyback infusion mode
6. RUN program or put on HOLD
7. Prime lines
8. Door latch
9. Clear programming
10. Tubing guide

MACRO and MICRO Infusion Setup

1. Plug pump into AC wall outlet or use on battery power.
2. Open front door and insert IV line into guide groove so that fluids flow from left to right. Close door.
3. Hold **SELECT** down, then quickly press and release **ON/OFF.** Release **SELECT**. [SLCT] appears on display.
4. Press **SELECT** until MACRO or MICRO indicator turns green, then press **RUN/HOLD**.
5. If needed, press **SELECT** to scroll to desired pressure setting (HI, NOR, LO) then press **RUN/HOLD**.
6. RATE indicator turns green. Press ▲ or ▼ to scroll to desired rate of administration, then press **RUN/HOLD**.
7. You may now set either the total time of administration or the volume to be infused. Press **SELECT** to highlight either TIME or VTBI, then use ▲ or ▼ to scroll to desired setting. Press **RUN/HOLD**.
8. Press **PRIME** to fill IV lines with fluid.
9. Press **RUN/HOLD** to begin delivery.
10. To stop delivery, press **RUN/HOLD**.

Checking Settings

- To check time remaining, press **SELECT** once.
- To check remaining volume to be infused, press **SELECT** twice.
- To check volume that has been infused, press **SELECT** 3 times.

Adjusting Settings

- Press **RUN/HOLD**.
- Press **SELECT** until parameter to be adjusted is displayed.
- Use ▲ or ▼ to adjust setting.
- Press **RUN/HOLD** to resume operation.

Other Infusion Options

10–P: Sets up to 10 different rate/volume combinations to be given to a single patient over time.

PIGY: Allows two separate rate/volume settings to be programmed to accommodate a main fluid bag plus a piggybacked supplemental bag.

TPN: Programs ramp-up, ramp-down and total delivery time for use with Total Parenteral Nutrition.

See Heska Vet/IV™ 2.2 Infusion Pump User Manual for details.

©2022 Heska Corporation. All Rights Reserved. Baxter is a trademark of Baxter International, Inc. Hospira is a registered trademark of Hospira, Inc. Zoetis is a registered trademark of Zoetis Products, LLC. Covetrus is a registered trademark of VFC Central, LLC. HESKA is a registered trademark and Vet/IV is a trademark of Heska Corporation in the U.S. and other countries. US22MD0101 43647-001